


# BAINLINE - H

Metallography  
Heavy-Duty Belt grinder

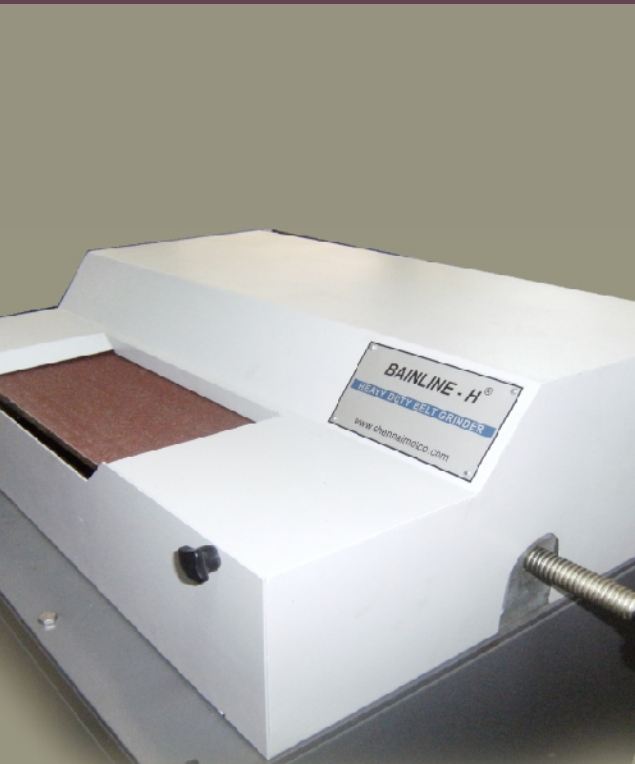

For flat, coarse grinding prior to fine grinding. The specimens to be ground in 60,80 or 120 grit after cutting operation or moulding to remove surface imperfections. Useful step to remove light metallurgical damage in the specimen surface induced during cut-off operation..

- For Pre-Grinding of Metallurgical Specimen & Large metal section
- Rigid Floor Model construction
- Variable Speed (Mechanical)
- Less Vibration Design
- 1 HP 3 phase motor
- 150mm x 1220mm Endless Abrasive Belt
- Constant Tension Spring Loaded Mechanical system
- Aluminium Pulleys with Belt Drive Mechanism
- Simplified Belt release mechanism
- Easy removable Dust cover system
- Size: 700x600x1120mm
- On-Off Starter (Optional)



**Bainline BG (Table Top) Single/Twin,  
and Heavy Duty Spectro Grinder available**



## **Chennai Metco Pvt Ltd**

SP 100A, Ambattur Industrial Estate, Chennai - 600 058, India.  
Ph: 044-2688 0737, 044-2688 0005 Fax: (044) 2688 0850,  
e-mail: [mktg@chennaietco.com](mailto:mktg@chennaietco.com), web : [www.chennaietco.com](http://www.chennaietco.com)