


Bainmount - Hydraulic

Hydraulic Hot Mounting Press

vBainmount - Hydraulic Metallography specimen mounting press is suitable for the small metallography specimen. After forming, it will be very easy to process the grinding, polishing of the specimen and it is also convenient to process measuring of metallographic structure under metallographic microscope.

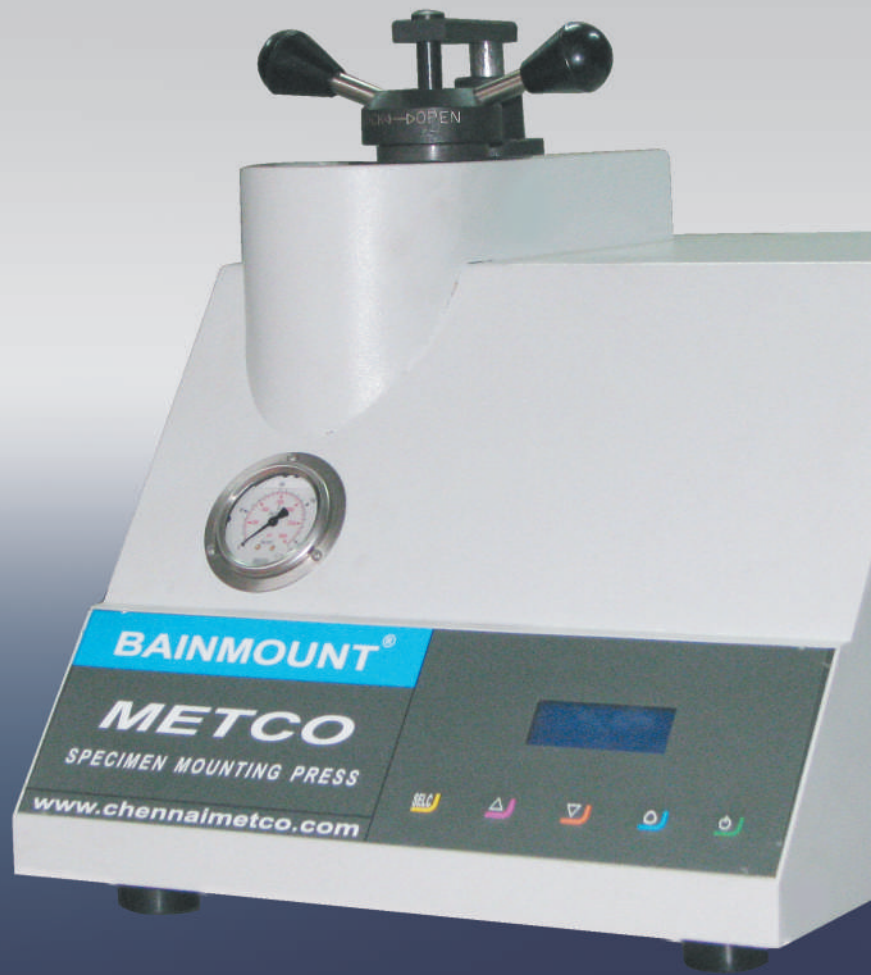

High Speed Heater, Direct Reading pressure gauge, Timer with Buzzer at end of cycle.

Chennai Metco



- Corrosion resistant steel Cabinet with built in Hydraulic system
- Individual pump and Cylinder
- Manual operation system
- Digital Temperature Indicator and digital timer
- Stainless steel body Pressure Gauge system
- Quick Ram return to speed up reloading
- 1500 Watt mould Heater
- Dovetail slot for easy changeability of Mould Assemblies
- Water Cooling System (manual)
- Timer-Buzzer indication for heating and cooling cycles
- Wide range of sample sizes can be accommodated
- Two moulds can be simultaneously prepared using Dual Mould Spacer
- Provision to accommodate 1 ½ " & 2" Mould Sizes
- Size: 490 x 620 x 590 mm

Optional Accessories

- Mould Assembly 1 ½ "
- Mould Assembly 2"
- Water Cooling (Automatic)

Pneumatic also available

Chennai Metco Pvt Ltd

SP 100A, Ambattur Industrial Estate, Chennai - 600 058, India.
Ph: 044-2688 0737, 044-2688 0005 Fax: 044 2688 0850,
e-mail: mktg@chennaietco.com, web : www.chennaietco.com