


Bainmount - H Auto

Automatic Hydraulic Hot Mounting Press

- Fully Auto Electro Hydraulic Mounting Press
- Pressure, Time and Temperature can be set to suit the cycle
- Fully automatic with auto cooling adjustable
- Useful for all types of mould materials
Bakelite, Transparent, Epoxy etc.

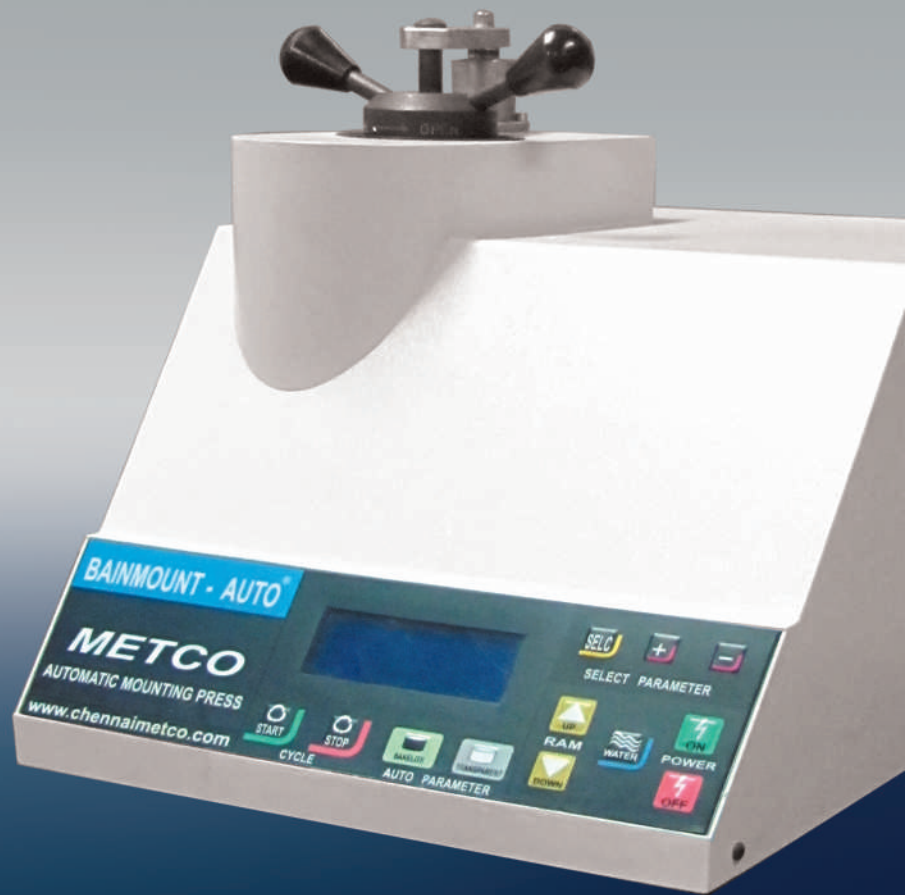

HIGH SPEED HEATER, DIRECT READING PRESSURE GAUGE, TIMER WITH BUZZER AT END OF CYCLE.

Chennai Metco

- Corrosion resistant steel Cabinet with built in Hydraulic system
- Individual power pack and Cylinder
- Hydraulic Fully Automatic system
- Digital Temperature Indicator and digital timer with LCD Display System
- Pressure transmitter (to change pressure)
- Quick Ram return to speed up reloading
- 1600 Watt mould Heater
- Dovetail slot for easy change of Mould Assemblies
- Water Cooling System (Automatic)
- Timer-Buzzer indication for heating and cooling cycles
- Wide range of sample sizes can be accommodated
- Two moulds can be simultaneously prepared using Dual Mould Spacer
- Provision to accommodate 1 ½ " & 2" Mould Sizes
- Size: L:600 X B:440 mm X H: 460 mm
- 1¼" comes with the machine as a standard
- Quick Chamber Open/close.

Optional Accessories

- Mould Assembly 1 ½ "
- Mould Assembly 2"
- Recirculating water cooling system of 30 Ltr capacity

Hydraulic Manual also available

Chennai Metco Pvt Ltd

SP 100A, Ambattur Industrial Estate, Chennai - 600 058, India.
Ph: 044-2688 0737, 044-2688 0005 Fax: (044) 2688 0850,
e-mail: mktg@chennaietco.com, web : www.chennaietco.com