


BAINCUT

XXL

AUTO

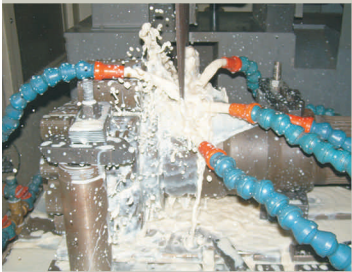

Automatic bigger size fully loaded cutting machine with top closed design. Voluminous internal space to keep big components.

www.chennai.netco.com

- Sturdy floor model cutter for sectioning Metallography samples
- Rust Proof Steel Construction
- Cutting Capacity upto 180 mm dia or W * H
- 25 HP, 3 Phase Motor
- AC Drive Speed Control
- Spindle speed Upto 100 - 2880 RPM
- Heavy Duty Steel construction
- Cut-off wheel dia 20"
- PLC Encoder With Control for 'Z' Movement Depth of cutting
- Cutting action automatic – ' Y ' Movement by Servo System
- Cutting action automatic-'Z' Movement by Servo system
- Splash proof, corrosion resistant Doors with see-through System
- Cooling by Multiple high flow water jets (5 numbers) to provide optimum cooling
- Coolant Tank with 200 ltrs capacity
- 1/2 HP Recirculation pumps (2 numbers)
- Two large drains for draining the coolant
- Flash cleaning attachment with hook for mounting
- Fluorescent light in the working area provides safe & clear illumination
- Start button with key for safety



XXL - PANEL



1. MMI
2. Joystick for Y Movement - Manual
3. Joystick for Z Movement - Manual
4. Y Movement Cycle Start
5. Z Movement Cycle Start
6. Jog Cycle Start
7. Lamp On/Off
8. Emergency Stop

Baincut — XXL Manual also available

Chennai Metco Pvt Ltd

SP 100A, Ambattur Industrial Estate, Chennai - 600 058, India.
Ph: 044-2688 0737, 044-2688 0005 Fax: (044) 2688 0850,
e-mail: mktg@chennaietco.com, web : www.chennaietco.com